

### **AYURVEDIC MASSAGE BED ( DHARAPATHI / DHRONI )**

1. FIBER - EXTRA WIDE
2. FIBER - STANDARD WIDTH
3. FIBER - COMPACT HEAVY DUTY
4. WOOD - TRADITIONAL MODEL
5. WOOD - SHORT BOWL MODEL
6. WOOD - RECTANGULAR MODEL

### **SUPPORT STANDS**

7. STEEL - HI-LO TILT TABLE
8. STEEL - TILT TABLE STAND
9. STEEL - STAND
10. FIBER - TULSI THARA STAND
11. WOOD - PLAIN STAND
12. WOOD - DECORATIVE STAND

### **AYURVEDIC MASSAGE TABLE ( SHIRODHARA TABLE )**

13. FIBER - STRETCHER TABLE
14. FIBER - ECONOMY TABLE
15. FIBER - LIGHT WEIGHT TABLE
16. WOOD - ABHYANGA TABLE
17. WOOD - DECORATIVE TABLE
18. WOOD - ALL AROUND TABLE
19. WOOD - OPEN SHELF TABLE
20. WOOD - DECORATIVE CABINET TABLE

### **PREMIUM DESIGNER MODELS**

### **MASSAGE ACCESSORIES**

21. STEEL - STEP
22. WOOD - STEP
23. WOODEN MASSAGE STOOL
24. SOFT BED
25. PILLOWS FOR HEAD & KNEES

### **DISPOSABLE THERAPY ACCESSORIES**

26. DISPOSABLE TREATMENT CLOTHING
27. OTHER DISPOSABLES

### **SHIRODHARA EQUIPMENTS**

28. WOOD - ORDINARY STAND
29. WOOD - DECORATIVE 2 LEG STAND
30. WOOD - DECORATIVE 4 LEG STAND
31. STEEL - FULLY ADJUSTABLE STAND
32. OIL HEATING SYSTEM
33. SHIRODHARA MACHINE FULLY AUTOMATIC
34. SHIRODHARA MACHINE SEMI AUTOMATIC
35. SHIRO DHARA VESSELS

### **STEAM BATH CHAMBER - SITTING TYPE**

36. WOOD - STANDARD MODEL
37. WOOD - ECONOMY MODEL
38. WOOD - DECORATIVE MODEL
39. FIBER - BOX MODEL
40. FIBER - CURVED MODEL

### **STEAM BATH CHAMBER - LYING TYPE**

41. WOOD - BOAT SHAPE ( BUTTERFLY COVER )
42. FIBER CHAMBER WITH STEEL FRAME
43. FIBER CHAMBER WITH WOOD FRAME
44. WOOD - HEXAGONAL MODEL
45. WOOD - TRAPEZIUM MODEL
46. WOOD - EXOTIC DESIGN

### **STEAM GENERATORS / BOILERS**

47. ELECTRIC BOILER
48. GAS BOILER
49. SEMI AUTOMATIC ELECTRIC BOILER
50. GAS & SEMI AUTO ELECTRIC BOILER (DUAL)
51. FULLY AUTOMATIC ELECTRIC BOILER
52. FULLY AUTOMATIC WALL MOUNTING TYPE
53. NADI SWEDAN YANTHRAM - HAND HELD
54. NADI SWEDAN YANTHRAM

### **LYING STREAM BATH WITH DHRONI**

55. MULTI PURPOSE PANCHAKARMA STATION

### **KSHARA SUTRA**

56. KSHARA SUTRA DRIER CUM STERILIZER

### **PIZHICHIL / SARVANGHA DHARA**

57. AUTO PIZHICHIL / SARVANGHDHARA MACHINE

### **AYURVEDIC BATH TUBS (AWAGAHA EQUIPMENTS)**

58. TUB WITH SUPPORT STAND
59. TUB WITH HEATING & CIRCULATION SYSTEM
60. SPINE BATH
61. HIP BATH ( KATI SNAN)

### **PANCHAKARMA EQUIPMENTS ( PEETH )**

62. VAMAN PEETH
63. VIRECHAN PEETH ( VIRECHAN CHAIR )
64. NASYA PEETH ( NASYA CHAIR )
65. ARSH BASHPA SWEDA YANTHRAM  
( RECTUM STEAM CHAIR )
66. RAKTHA MOKSHAN CHAIR ( LEECH THERAPY )

### **AYURVEDIC SAUNA ( THAAP SWEDAN YANTHRA )**

67. SITTING MODEL
68. LYING MODEL

### **VASTI YANTHRAMS ( BASTI )**

69. VASTI YANTRAM
70. MATHRA VASTI / THAILA VASTI YANTRAM
71. NIRUHA BASTI - DISPOSABLE SET
72. UTTAR BASTI YANTRA
73. BAHYA BASTI YANTRA

### **OTHER THERAPY ACCESSORIES**

74. PANCHAKARMA ACCESSORIES
75. SHALAKA ( PANCHALOHA )
76. SHIROVASTI CAP
77. FACIAL STEAMER / VAPORIZER
78. IR HEATER WITH FLAT URULI

### **SPA FURNITURE**

79. SPA BED NEW
80. MASSAGE TREATMENT TABLE
81. KANSA VATKI BOWL & WAND
82. FOOT BATH BOWL
83. WOODEN MULTI PURPOSE BENCH \ OIL VESSEL
84. WOODEN TROLLEY

### **THERAPEUTIC PREPARATIONS**

85. AYURVEDIC MEDICINES
86. AYURVEDIC / SPA OILS
87. HERBAL PREPARATIONS

## 1 FIBER - EXTRA WIDE

Export Model



- Heavy duty extra wide fiberglass dharapathi.
- Made with superior NPG- ISO quality FRP resin.
- Has double reinforcement to carry any weight.
- Filled with polyurethane foam (PUF) to resist impacts during handling/ transportation / shipping.
- Standard model is available in plain colours.
- Export & customized model available in teakwood / rosewood / jungle wood design surface finish.

Can be placed on Steel /Wooden/Fiber Support Stands

cm : 290 x 90 x 13  
inch : 114 x 36 x 5

## 2 FIBER - STANDARD WIDTH

Standard Model



- Heavy duty fiberglass (FRP/GRP) dharapathi of standard size.
- Made of superior NPG- ISO quality FRP resin.
- Assured life-long performance without any erosion, corrosion or need for polishing.
- Available in wood design surface finish as well as plain colours.

Can be placed on Steel / Wooden Support Stands

### 8 Feet Model

cm : 247 x 78 x 11  
inch : 97 x 31 x 4

### 9.5 Feet Model

cm : 290 x 79 x 13  
inch : 114 x 31 x 5

## 3 FIBER - COMPACT HEAVY DUTY

Export Model



- Heavy duty FRP dharapathi of compact size.
- Made of superior NPG- ISO quality FRP resin.
- Though overall length is shorter , the treatment area is same as that of an 8 feet dharapathi.
- Suitable for small treatment rooms.

Can be placed on Steel/Wooden/Fiber Support Stands

cm : 220 x 79 x 13  
inch : 87 x 31 x 5

## 4 WOOD - TRADITIONAL MODEL

Export Model



- Heavy duty wooden dharapathi.
- Made of country venga / other medicinal wood.
- Kiln & oil seasoned with steam / herbal oil for enhanced strength, durability & crack-resistance.
- Exquisitely hand crafted design with decorative brass work.

Can be placed on Steel/Wooden/Fiber Support Stands

### Normal Width

cm : 280 x 76  
inch : 110 x 30

### 3 Feet Width

cm : 280 x 90  
inch : 110 x 35

## 5 WOOD - SHORT BOWL MODEL

Export Model



- Has all the features of 4 Wood -Traditional Model.
- Head portion of the dharapathi has a head bowl.
- This increases the crack resistance & provides more space to therapist for head massage.

Can be placed on Steel / Wooden Support Stands

### Normal Width

cm : 280 x 76  
inch : 110 x 30

### 3 Feet Width

cm : 280 x 90  
inch : 110 x 35

## 6 WOOD - RECTANGULAR MODEL



- Extra wide dharapathi with a space saving design.
- Made of medicinal wood.
- Provided with a slope towards leg –end.
- Allows even a 6.5 feet tall person to have comfortable abhyanga / shirodhara treatment with only 8 feet long table.

Can be placed on Steel / Wooden Support Stands

cm : 240 x 90  
inch : 94 x 35

# Support Stand for Massage Bed / Dharapathi / Dhroni

## 7 STEEL - HI-LO TILT TABLE

With Tilt Model



- Simple & efficient mechanism to raise & lower the table to near ground level.
- Allow oiled patients who are fat / have large body frame to climb up & get down from the dharapathi safely & effortlessly without the help of an assistant.
- Tilt model enables tilting of the head-end & leg-end for shirodhara / anjana / pizhichil / sarvangadhara procedures.
- All types of dharapathi can be placed on this table.

Prices quoted are for the Stand only.  
Dharapathi sold seperately

## 8 STEEL - TILT TABLE STAND




- Innovative design to tilt the table top.
- Rotating handle wheel clockwise lowers the head-end for shirodhara procedures.
- Rotating handle wheel anti-clockwise lowers the leg-end for pizhichil / sarvangadhara procedures.
- Sturdy construction with 6 legs and two tie frames of powder coated steel.

# Support Stand for Massage Bed / Dharapathi / Dhroni

## 9 STEEL - STAND

Powder Coated Steel Stand

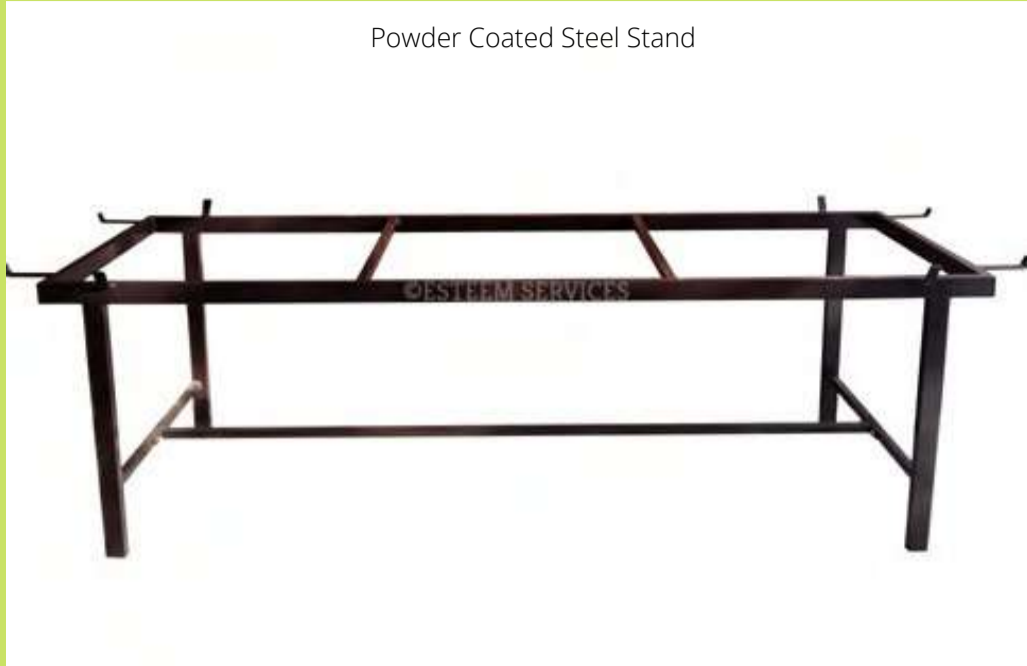

- Less heavier than **8** Steel - Tilt Table Stand.
- Made of steel pipes that are electrostatically powder coated for corrosion protection.
- Can also be made with stainless steel material.
- Collapsible design is provided to reduce freight cost by dismantling the legs.

**8 Feet Model**

**9.5 Feet Model**

## 10 FIBER - TULSI THARA STAND

Standard Model



- Attractively shaped as 'tulasi thara'.
- Ballasted for added stability.
- Allows maximum leg room around the dharapathi, compared to other stands & tables.
- Available in plain and wood designs to match Dharapathi

Prices quoted are for the Stand only.  
Dharapathi sold separately

## 11 WOOD - PLAIN STAND

Standard Quality Model



- Excellently polished strong frame with designed legs.
- Made of medicinal wood / hard wood.
- Provided with inter-tied legs & brass decorations.
- All types of dharapathi can be placed on this table

## 12 WOOD - DECORATIVE STAND



- Highly decorative wooden stand with hand carved floral designs on a 6 inch size wooden frame.
- Made out of country vengal/ other medicinal wood.
- Provided with brass work and 4.5 inch size inter tied legs.
- Brass flowers are provided on the floral designs.

## 13 FIBER - STRETCHER TABLE ( PORTABLE )

Standard Model



- Strong and lightweight for multipurpose use.
- Made of high strength fiberglass material.
- Placed on it's own special (collapsible) table for portable use.
- Invalids & paraplegics can undergo abhyangam / shirodhara/ pizhichil treatment in their own rooms.
- These patients can also be carried on this dharapathi for other treatments, bath etc.

cm : 274 x 74 x 8  
inch : 108 x 29 x 3

## 14 FIBER - ECONOMY TABLE ( PORTABLE )

Wood Design Model



- Economical dharapathi intended for new startup clinics.
- Placed on it's own special (collapsible) table for portable use.
- The legs can be removed & re-assembled easily for carrying the table to patient's place.
- Can be used for shirodhara / pizhichil / abhayanga procedures.

7 Feet Model

cm : 207 x 91 x 6  
inch : 81 x 36 x 2

8 Feet Model

cm : 235 x 81 x 6  
inch : 93 x 32 x 2

## 15 FIBER - LIGHT WEIGHT TABLE



- Easily portable lightweight dharapathi.
- Made from high strength FRP raw materials.
- Has traditional pathy's shape of head & drain
- The dharapathi has to be used with its own steel stand which is made of rectangular powder coated steel pipes.
- Legs can be removed easily.
- Available in many colour options.

cm : 274 x 74 x 8  
inch : 108 x 29 x 3

## 16 WOOD - ABHYANGA TABLE



- Strong heavy duty frame.
- Made of medicinal wood.
- Top surface has a slant towards the leg end to facilitate easy draining of massage oils.
- Can be used for abhyangam / pizhichil / vasti procedures.

cm : 195 x 75 x 75  
inch : 78 x 30 x 30

## 17 WOOD - DECORATIVE TABLE



- Export wooden Dhroni is placed on a highly decorative wooden stand with 6 inch framework
- Dhroni is made of medicinal woods and seasoned.
- Brass flowers are provided on the hand-carved floral designs.

cm : 270 x 75 x 75  
inch : 108 x 30 x 30

cm : 270 x 90 x 75  
inch : 108 x 36 x 30

## 18 WOOD - ALL AROUND TABLE



- Export quality wooden dharapathi.
- Highly decorative wooden table.
- Both are made of medicinal wood.
- The table has a 6 feet size frame with beautiful hand-carved floral designs encircling the dhroni, brass flowers and decorative work.
- Provided with brass decorations and inter - tied legs.

## 19 WOOD - OPEN SHELF TABLE



- Compact massage therapy table (8 feet long) with bowl-less dharapathi.
- The straight ends of dharapathi allow therapists a closer access for massage of head & legs.
- The dharapathi is polished & stained to match the colour of table.
- Support table has open shelves for storage of medicines & accessories.
- Can also be used for abhyangam/ shirodhara / pizhichil treatments.

Massage table accessories like soft bed, steps, panchakarma utensils (bronze) etc. are also shown in this brochure.

## 20 WOOD - DECORATIVE CABINET TABLE



- Export quality wooden dhroni placed on highly decorative cabinet table.
- Holders are provided for oil vessels on both ends, to collect & store used oils.
- Hand carved designs on 4 decorative doors on either side of the table.
- Lion leg shaped hand carved legs are also provided.
- Excellently finished with high quality sealers & poly urethane (PU) topcoat.

Massage table accessories like soft bed, steps, panchakarma utensils (bronze) etc. are also shown in this brochure.

# Premium Designer Models

# prices available on request

**A** WOODEN DHRONI ON OPEN SHELF TABLE ( Natural Anjali )



**B** WOODEN DHRONI ON OPEN SHELF TABLE (Walnut Finish)



**C** WOODEN DHRONI ON OPEN SHELF TABLE WITH DRAWER



**D** RECTANGULAR DHRONI ON SHELF TABLE



**E** WOODEN DHRONI ON TABLE WITH OIL VESSEL DRAWERS



**F** DECORATIVE HALF CABINET TABLE



**G** WOODEN DHRONI ON BANYAN TREE BASE



**H** WOODEN DHRONI ON DECORATIVE TABLE



# Premium Designer Models

# prices available on request

**I** FIBER DHRONI ON DECORATIVE ALL AROUND TABLE



**J** WOODEN MASSAGE TABLE WITH FACE HOLE & LID



**K** FIBER EXTRA WIDE DHRONI ( Cedar Wood Finish )



**L** FIBER DHRONI ( Jungle Wood Finish )



**M** SHORT BOWL WOODEN SHIRODHARA TABLE



**N** SHORT LENGTH DHRONI WITH STORAGE PLANK



**O** SHIRODHARA TABLE WITH ALL AROUND CARVINGS



**P** STEAM CHAMBER LYING BOX TYPE





# Accessories for Ayurvedic Massage Bed / Dharapathi

## 21 STEEL - STEP



### Mild Steel Step:

- Made from high quality CR pipes.
- Electrostatically powder coated for long life.
- Single step and double stair are available

### Stainless Steel Step:

- Made from superior quality stainless steel pipes

## 22 WOOD - STEP



- Made from seasoned wood.
- Hand carved design work is made in Export quality step.

## 23 WOODEN MASSAGE STOOL



- Wooden massage stool made of good quality wood.
- Aesthetically carved legs and excellent polish are the plus points.
- Essential for head oiling preceding Abhyangam.

# Accessories for Ayurvedic Massage Bed / Dharapathi

## 24 SOFT BED



- The soft bed relieves distress of the patient lying on a hard dharapathi.
- Made of sheathed superior quality foams enclosed in zipped artificial leather covers to resist oil soaking.
- Top part is shaped to match the curvature of the shoulder hump of the dharapathi.
- Allows replacement of the inner foam when required.

## 25 PILLOWS FOR HEAD & KNEES



- Two Knee pillows and a head/face pillow reduce discomfort at the knee joints and jaws during massage in prone position.
- Specially made to resist oil soaking using superior quality foams enclosed in zipped artificial leather covers.
- This allows replacement of the inner foam when required.

# Disposable Therapy Accessories

## 26 DISPOSABLE TREATMENT CLOTHING



Economical disposable items made of skin friendly nonwoven material. Shown here in clockwise direction as below :

- Apron (Disposable)
- Skirt (Disposable)
- T Band (Ladies & Gents)
- Bra Elastic (Medium & Large)
- Panties Elastic (Small, Medium, Big).

Detailed Price list available on request.

## 27 OTHER DISPOSABLES



Economical disposable items made of skin friendly nonwoven material. Shown here in clockwise direction as below :

- Patient's gown
- Dhroni (bed) sheet
- Pillow covers (Small, Medium, Large).
- Shirodhara band
- Shower cap
- Facial band
- Head band
- Eye cover
- Hand gloves
- Face masks.

Detailed Price list available on request.

# Shirodhara Stands

## 28 WOOD - ORDINARY STAND



- Simple shirodhara stand economical model.
- Made of hardwood with 4 legs.

Prices quoted are for the Stand only.  
Dhara Vessel sold separately.

cm : 170 x 55  
inch : 67 x 22

## 29 WOOD - DECORATIVE 4 LEG STAND



- Two legged stand with adjustable hook & wheels.
- Made from medicinal wood.
- Legs go under the table saving floor space.
- Collapsible design to reduce shipping costs.
- Provided with beautiful hand carvings.

Prices quoted are for the Stand only.  
Dhara Vessel sold separately.

cm : 170 x 55  
inch : 67 x 22

## 30 WOOD - DECORATIVE 2 LEG STAND



- Decorative wooden shirodhara stand.
- Hand-carved out of medicinal wood.
- Collapsible design to reduce shipping costs.
- Provided with beautiful brass decorations.

Prices quoted are for the Stand only.  
Dhara Vessel sold separately.

cm : 186 x 55  
inch : 73 x 22

## 31 STEEL - FULLY ADJUSTABLE STAND



- Sturdy steel shirodhara stand.
- Made in heavy duty stainless steel / rust-resistant fully powder coated mild steel pipes.
- Ballasted for perfect balance & stability.
- The shirodhara vessel can be positioned anywhere in both horizontal & vertical directions.

Prices quoted are for the Stand only.  
Dhara Vessel sold separately.

cm : 180 x 55  
inch : 70 x 22

# Shirodhara Equipments

## 32 OIL HEATING SYSTEM (OHS)



- Consists of a traditional hand made white bronze vessel with an attached heating coil & an energy controller.
- Provides a steady flow of shirodhara oil at constant heat.
- Consumes very low electrical power.
- Can be used with any steel / wooden shirodhara stand.

Prices quoted are for the OHS only.  
Stand sold separately.

## 33 SHIRODHARA MACHINE FULLY AUTOMATIC



- Portable unit with a steel shirodhara stand & a shirodhara machine positioned at its base.
- Instead of pumping oil directly on forehead (which increases pressure of shirodhara liquid) , it's pumped to moving shirodhara vessel & falls on the forehead .
- The machine filters used oil & warms it to constant temperature before silently pumping it to the oscillating vessel positioned over patient's head.
- The swinging width & speed can be adjusted.
- The vessel can also be moved manually without power.
- Provided with digital temperature controller & a machine body heater for cold climates.

## 34 SHIRODHARA MACHINE SEMI AUTOMATIC



- This is similar to **29** Shirodhara Machine Fully Automatic model.
- However the shirodhara vessel needs to be moved manually.

# Shirodhara System

## 35 SHIRO DHARA VESSELS

All these vessels come with Brass Chains, Rings, hook & imported ball valve control or cotton wick.



### FIBERGLASS MODELS ( FRP / GRP )

- FRP Vessels are very strong, lightweight & easy to clean.
- These are economical vessels that last a lifetime.
- The wood design vessels are made in teakwood / rosewood / junglewood shades to match the dhroni or steam chamber.
- Provided with stainless steel chains in plain design model & brass chains in wood design model.



### COPPER MODELS

- These are made of high purity copper sheets.
- The vessels are lightweight and easy to handle.
- The copper pot model has stainless steel chains and the vessel model comes with brass chains.



### WHITE BRONZE MODELS

- The white bronze metal is traditionally used to make vessels & other apparatuses in kerala ayurveda practices.
- These are considered superior to any other metals in commercial use.



### OTHER MODELS

#### BRONZE VESSELS

- They retain heat and are more resistant to the actions of herbs & water compared to the brass.
- Brass not used for utensils as per traditional ayurveda practices.

#### STAINLESS STEEL VESSELS

- These are made with good quality stainless steel.
- Also provided with stainless steel hanging chains.

## 36 WOOD - STANDARD MODEL



Width of the door of the treatment room should be at least 84 cm to suit the width of chamber.

- Heavy duty 4 doors wooden steam chamber.
- Made from medicinal woods / anjali.
- Provided with copper tube steam release system.
- Standard model has polymer inner lining.
- Export model has fiber inner shell & molded non-slip flooring inside the wooden chamber.
- This protect the wood from steam & sweat & maintains a hygienic surface without any crevices where bacteria may lodge.
- Adjustable stool & brass fittings is also provided.
- Digital temperature controller is optional ( ₹ 4800 extra ).

cm : 90 x 82 x 137  
inch : 35 x 32 x 54

## 37 WOOD - ECONOMY MODEL



Width of the door of the treatment room should be at least 76 cm to suit the width of chamber.

- Heavy duty 2 doors wooden steam chamber.
- Made of wood with water proof interior surface.
- Heat resistant steam dispenser is provided.
- Adjustable stool with stainless steel top is also provided.

cm : 110 x 74 x 136  
inch : 43 x 29 x 54

## 38 WOOD - DECORATIVE MODEL



Width of the door of the treatment room should be at least 84 cm / 92 cm to suit the width of chamber.

- Decorative export quality 2 door steam chamber.
- Made from medicinal woods / anjali.
- Provided with fiber (FRP) inner lining & molded non-slip flooring inside the wooden chamber
- This protect the wood from steam & sweat and maintains a hygienic surface without any crevices where bacteria may lodge.
- Doors have aesthetically hand carved traditional designs.
- Adjustable wooden seat is provided.
- Digital temperature controller is optional ( ₹ 4800 extra ).

cm : 90 x 82 x 137  
inch : 35 x 32 x 54

cm : 90x90 x 137  
inch : 35 x 35 x 54

## 39 FIBER - BOX MODEL



Width of the door of the treatment room should be at least 68 cm / 78 cm to suit the width of chamber.

- Fiber (FRP) heavy duty steam chamber.
- High temperature FRP resin (HDT > 100 degree C) assures life long performance without need for maintenance.
- Provided with heat resistant FRP claddings for copper pipes & moulded non-slip floor.
- Provided with copper tube steam release system.
- Adjustable seat & steel support stand is provided.
- Digital temperature controller is optional ( ₹ 3900 extra ).

cm : 93 x 66 x 116 / cm : 93 x 76 x 116  
inch : 37 x 26 x 46 / inch : 37 x 30 x 46

## 40 FIBER - CURVED MODEL



Width of the door of the treatment room should be at least 82 cm to suit the width of chamber.

- Heavy duty fiber steam chamber with solid box type door.
- Chamber surface is finished in junglewood / rosewood / teakwood designs or plain colours.
- High temperature FRP resin (HDT > 100 degree C) assures life long performance without need for maintenance.
- Provided with copper tube steam release system.
- Heat resistant FRP claddings for copper pipes, moulded non-slip floor & PUF insulated door are provided.
- Digital temperature controller is optional ( ₹ 3900 extra ).

cm : 92 x 80 x 124  
inch : 36 x 31 x 49

## 41 WOOD - BOAT SHAPE ( BUTTERFLY DESIGN )



Width of the door of the treatment room should be at least 80 cm to suit the width of chamber.

- Ideal steamer for Spa.
- Space saving design in attractive boat shape.
- Bed is made of medicinal wood.
- Window panes are made of acrylic glass.
- FRP steam condensate tray prevents steam from entering the lower wooden chamber. The steam circulates within the tray & upper chamber.
- This protects the wooden surface & increases life of wood.
- Provided with stainless steel steam distributors.
- Digital temperature controller is optional ( ₹ 4800 extra ).

cm : 225 x 78 x 100  
inch : 89 x 31 x 39

## 42 FIBER CHAMBER WITH STEEL FRAME

Powder Coated Steel Model



Width of the door of the treatment room should be at least 80 cm to suit the width of chamber.

- Bed is made of teakwood / medicinal wood.
- Provided with powder coated steel frame & 2 piece bed.
- Chamber is made in translucent / opaque heat resistant fiberglass.
- Patient can be transferred to the bed from the sides by lifting-up the FRP cabin.
- Provided with stainless steel tube steam release system.
- FRP condensate tray & foldable headrest are also provided.
- Digital temperature controller is optional ( ₹ 3900 extra ).

cm : 210 x 78 x 90  
inch : 83 x 31 x 35

## 43 FIBER CHAMBER WITH WOOD FRAME

Wooden Base Model



Width of the door of the treatment room should be at least 86 cm to suit the width of chamber.

- Base frame & bed are made of teakwood / medicinal wood.
- Chamber is made of translucent / opaque heat resistant fiberglass with stainless steel framework.
- Patient can be transferred to the bed from the sides by lifting-up the FRP cabin.
- Provided with stainless steel tube steam release system.
- FRP condensate tray & foldable headrest are also provided.
- This model is maintenance free & easy to clean.
- Digital temperature controller is optional ( ₹ 3900 extra ).

cm : 210 x 84 x 90  
inch : 83 x 33 x 35

## 44 WOOD - HEXAGONAL MODEL

Hardwood Frame Model



Width of the door of the treatment room should be at least 86 cm to suit the width of chamber.

- Bed is made of teakwood placed on a wooden stand made of medicinal wood.
- Chamber is made of wooden frame & translucent, wrinkled FRP window panes made with heat resistant NPG fiberglass.

cm : 210 x 84 x 88  
inch : 83 x 33 x 35

## 45 WOOD - TRAPEZIUM MODEL

Standard Model



Width of the door of the treatment room should be at least ? to suit the width of chamber.

- Economical version of 40 Wood - Hexagon Model.
- Similar to Wood - Hexagon Model in construction.
- Has a simpler chamber shape with hylam (opaque) or translucent FRP windows.

cm : 190 x 80 x 90  
inch : 75 x 32 x 36

## 46 WOOD - EXOTIC DESIGN

Lion Model



Elephant Model



Width of the door of the treatment room should be at least 92 cm to suit the width of chamber.

- Exotic whole body wooden steam bath.
- Considered an important ayurvedic spa equipment.
- Steam bath has a royal look with majestic statues of Lion / Elephant supporting the head and leg ends.
- Top cabin is made as 4 sections for easy opening by hand.
- The fiber shell provided inside the top chamber & FRP steam condensate tray in the lower chamber protect the wooden parts from steam & sweat.
- This helps to maintain a hygienic surface without any crevices where bacteria may lodge.
- Has wooden cabinets under the bed to store the Steam Boiler and all treatment accessories, medicines etc.

cm : 280 x 90 x 130  
inch : 110 x 35 x 51



## 47 ELECTRIC BOILER



- Electric boiler of 2 KW rating.
- Made of superior SS 304 stainless steel metal.
- High quality ISI rated heating coil of 2 KW capacity is used.
- Pressure guage, safety valve, weight valve, ball valve & steam outlet are also provided.
- Perfume basket for soaking herbs enable the generation of herbal steam.
- Water level check valve helps to judge the water levels before starting and hence ensure safety against damages to the heater.
- Provided with a powder coated stand.
- SS Perfume box, droplet remover, ball valve, pressure gauge and safety valve are provided.
- Supplied with powder coated steel stand.

cm : 36 x 36 x 60  
inch : 14 x 14 x 24

## 48 GAS BOILER



- This is a Stainless Steel boiler (steam generator) for manual operation using LPG Gas for heating.
- Supplied with Gas Burner, Fittings, Stand, SS Perfume Box, safety valve, 'Waree' Pressure gage, Neoprene 'O' ring, cast steel Ball valve and a special Steam Hose with end fittings for leak proof joints.

cm : 36 x 36 x 72  
inch : 14 x 14 x 28

## 49 SEMI AUTOMATIC ELECTRIC BOILER



- This is a stainless steel (superior SS-304 grade- for use with Herbs) boiler of 2 KW capacity, with fully insulated double shell body (in Export quality)
- Electronic level sensing system (ELSS ) puts off the boiler & sounds alarm when water level reaches the heating coil, to avoid dry heating and damage of the element. When water is re-filled the boiler starts heating.
- Perfume box, control box, droplet remover( all of Stainless steel), trolley with wheels, ball valve steam control, Baumer pressure gauge, German thermostat, safety valve, Neoprene 'O' ring are provided.
- Standard quality comes with uninsulated body..

cm : 38 x 38 x 60  
inch : 15 x 15 x 24

## 50 GAS & SEMI AUTOMATIC ELECTRIC BOILER ( DUAL )



- This is a stainless steel (superior SS-304 grade- for use with Herbs) boiler of 2 KW capacity, with fully insulated double shell body (in Export quality)
- Electronic level sensing system (ELSS ) puts off the boiler & sounds alarm when water level reaches the heating coil, to avoid dry heating and damage of the element. When water is re-filled boiler starts heating.
- Perfume box, control box, droplet remover( all of Stainless steel), stand, ball valve steam control, Baumer pressure gauge, German thermostat, safety valve, Neoprene 'O' ring are provided.
- Standard quality comes with uninsulated body..
- Supplied with Gas Burner, Gas Fittings, Stand and special Steam Hose.

cm : 36 x 40 x 72  
inch : 14 x 16 x 28

## 51 FULLY AUTOMATIC ELECTRIC BOILER



- Fully automatic Electric Boiler (FAB) of 3 KW rating.
- Made of Stainless Steel (SS-304 grade) with an integral Alarm Warning System which switches off the Boiler if the water level inside goes down to prevent Heater damage.
- Can be connected to water line to allow automatic intake of water when required.
- A Perfume box is provided to insert herbs inside the boiler, to produce herbal steam bath.
- Fully insulated to avoid heat loss and save electric energy.
- All the fittings bolts & nuts are of SS or Brass to resist corrosion.
- Comes with Digital temp. controller (DTC) which is optional. A wheeled trolley is provided for mobility in Export quality.

cm : 62 x 42 x 40  
inch : 24 x 17 x 16

## 52 FULLY AUTOMATIC WALL MOUNTING TYPE



- Fully automatic Electric Boiler (FAB) of 3 KW rating.
- Made of Stainless Steel (SS-304 grade).
- Fully insulated to avoid heat loss and save electric energy.
- Integral alarm warning system switches off the boiler if the water level inside goes down to prevent heater damage.
- Can be connected to water line to allow automatic intake of water when required.
- A perfume box is provided to insert herbs inside the boiler, and produce a herbal steam bath.
- All the fittings bolts & nuts are of SS or Brass to resist corrosion.
- A wheeled trolley is provided for mobility.
- Digital temperature controller is optional at extra cost.

cm : 62 x 42 x 40  
inch : 24 x 17 x 16

## 53 NADI SWEDAN YANTHRAM- HAND HELD



- This swedana equipment is a portable (handheld) model for localized non- pressurized steaming of joints, back etc.
- A special steam hose with handle and water trap, SS Perfume Box, safety valve, ISI rated heating coil etc. are provided.
- This is made of Stainless Steel (SS304) to assure long life.

cm : 32 x 29 x 22  
inch : 13 x 11 x 9

## 54 NADI SWEDAN YANTHRAM



- This swedana equipment is a portable, insulated (in export model), electrical steam generator, used for localized non- pressurized steaming of joints, back etc.
- A special steam hose with handle and water trap, SS Perfume Box, safety valve, water level check valve, imported thermostat, ISI rated heating coil etc. are provided.
- This is made of Stainless Steel (SS304) to assure long life. A well balanced Stainless Steel trolley is provided in export range.

cm : 36 x 36 x 60  
inch : 14 x 14 x 24

## Lying Stream bath with Dhroni

### 55 MULTI PURPOSE PANCHAKARMA STATION



- Space & cost saving equipment serves that can be used as both Steam bath & Dhroni (ayurvedic massage bed).
- The body is made of epoxy coated steel. Lightweight Dhroni is made of high strength FRP material.
- Stainless steel electric boiler (SS-304 quality) supplies pressurized steam to the steam box.
- Refill water in the boiler after every 3 steam sessions.
- Wooden bed and shirodhara stand with Shirodhara vessel with ball valve control are provided.

cm : 220 x 82 x 80  
inch : 87 x 32 x 31

## Kshara Sutra

### 56 KSHARA SUTRA DRIER CUM STERILIZER ( WITH STORAGE )



Automatic Model

- This is an electric drier with ultraviolet sterilizer to make Kshara Sutra threads.
- Drying time for each of the eight stages can be pre-set in Automatic model.
- A storage cabinet with additional sterilizer for dried threads and an application tray with brush are provided in higher models.
- The stored threads are sterilized with each opening of the door.

cm : 46 x 36 x 61  
inch : 18 x 14 x 24

## Pizhichil / Sarvanghadhara

### 57 AUTOMATIC PIZHICHIL / SARVANGHADHARA MACHINE

Standard Model



- This is portable, compact unit and needs minimum space for operation.
- The bottom housing moves below the Dhroni.
- This machine enables Sarvanga Dhara, Ekanga Dhara, Shirodhara, Pradesika Dhara & Pizhichil to be done with one or two masseurs only.
- This has re-circulation system, filtering, oil heating system. provided with Dhara dispensers and adjustable limit switches for fixing travel length of Dhara vessel.
- Shirodhara vessel is to be moved manually in standard model .
- Automatic Shiro Dhara unit is provided in export model for an extra cost.
- Designed as per guidelines from classical texts.' dhara kalpam of sahasrayogam.

cm : 265 x 35 x 184  
inch : 104 x 14 x 72

## 58 TUB WITH SUPPORT STAND

Standard Model



- This is a very useful Bath Tub used in the treatment of ailments of Rheumatic nature, Neuralgia, skin diseases etc.
- Self-supporting with a steel frame, this can be placed anywhere in a hospital.
- Light weight yet very sturdy.
- Made of heat resistant NPG-ISO fiberglass.
- Support frame is of powder coated steel.

cm : 136 x 77 x 45  
inch : 54 x 30 x 18

cm : 168 x 77 x 55  
inch : 66 x 30 x 22

## 59 TUB WITH HEATING & CIRCULATION SYSTEM

Export Model



- This is a bath tub with facilities for heating, re-circulation and body -jet spray of hot water/ medicinal water.
- The tub is of heat resistant NPG fiberglass with handles and hand switch control.
- Digital controller assures correct temperature of water.
- Tub is clad with corrosion -proof Fiberglass body.

cm : 170 x 78 x 55  
inch : 67 x 31 x 22

## 60 SPINE BATH

Export Model with Groove



- Essential equipment in treatment of infertility and ailments of the back.
- This FRP spine bath is sturdy. Available in Plain, Rose Wood or Jungle Wood finish with spine groove.
- Provided with steel reinforcement at the periphery, in higher models.

cm : 101 x 63 x 22  
inch : 40 x 25 x 9

## 61 HIP BATH ( KATI SNAN )

Standard Model



Export Model with Stand



- This Sitting type tub offers a comfortable seating and is useful in treatment of ailments of the Pelvic region, sciatica Lumbago etc.
- Made of light weight FRP material and hence easily portable.
- Made of superior NPG-ISO fiberglass.

cm : 92 x 74 x 51  
inch : 36 x 29 x 20

## 62 VAMAN PEETH

Powder Coated Steel Model



- This is designed for executing Vamanam procedure comfortably.
- The chair is designed to be stable in the inclined sitting position.
- Made out of Stainless or Powder coated Steel (PC).
- Height of the basin is adjustable.
- Transparent measuring jar, trays for BP apparatus and tray for cup are provided.
- Fitted with wheels for mobility. Brakes are provided on wheels for arresting movement.
- Cushion seat and back is available at extra cost.

cm : 90 x 90 x 80  
inch : 35 x 35 x 31

## 63 VIRECHAN PEETH ( VIRECHAN CHAIR )

Powder Coated Steel Model



- This equipment facilitates comfortable Virechan (Purgation) procedure.
- The Stainless steel collection vessel enables easy and hygienic disposing of the excreta.
- Available with SS (Stainless steel) or PC (Powder coated) steel frame.
- Castor wheels enable mobility.
- The back is adjustable to any reclining position.
- Seat is made of wood design FRP.

cm : 75 x 57 x 81  
inch : 30 x 22 x 32

## 64 NASYA PEETH (NASYA CHAIR)

Powder Coated Steel Model



- This is a fully reclining chair with adjustable head rest, for Nasya & Anjana treatments as well.
- It can also be used as a lounge chair.
- The back rest and seat have teakwood/medicinal wood members.
- The frame can be of PC (powder coated) steel or stainless steel (SS).
- A new model with raisable leg support is also available for performing foot massage during Nasya Karma.

cm : 102 x 64 x 90  
inch : 40 x 25 x 35

## 65 ARSH BASHPA SWEDA YANTRAM ( RECTUM STEAM CHAIR )

Stainless Steel Model



- This is very useful equipment for treatments of Ano-rectal diseases like Piles, Fistula etc.
- This is provided with a padded back rest, arm rest & Wooden/FRP seat.
- Electric heater boils the water.
- A control box is provided below the arm rest with a heat control system and a regulator controlling the steaming period.

cm : 57 x 64 x 90  
inch : 22 x 25 x 35

## 66 RAKTHA MOKSHAN PEETH (LEECH THERAPY CHAIR )



- This is a comfortable tall chair for seating the patient with legs resting on the footrest.
- Seat height is adjustable to keep the legs straight
- Frame & legs are made from strong steel pipes..

cm : 47 x 44 x 92  
inch : 18 x 17 x 36

## Ayurvedic Sauna

### 67 SITTING MODEL



- This enables hot air sauna for single person.
- The wooden cabin comes with glazed window, Cabin light, Wooden Seat, Digital temperature indicator & controller, Floor mat etc. as standard accessories.
- The hot air generator is of powerful 3KW rating.
- A mechanical thermostat is provided as an additional safety feature.

cm : 75 x 75 x 180  
inch : 30 x 30 x 71

## Thaap Swedan Yanthra

### 68 LYING MODEL



- This equipment circulates dry, hot air around the person lying inside the chamber with head exposed outside.
- Digital temperature controller and a mechanical thermostat ensure safe, controlled heating.
- Bed is made of Teakwood or medicinal wood.
- Chamber cover is made of NPG fiberglass in export quality wood finish exterior.
- The structural part is made of power coated steel covered by FRP or similar insulating claddings.

cm : 212 x 85 x 90  
inch : 84 x 33 x 35

## 69 VASTI YANTRAM



1. Hand casted in the traditional way & shape by keralite artisans using authentic White bronze metal and has smooth, round tipped Nethram of 5 mm dia.
2. The one piece vasti yanthram is made of surgical plastic, for single person use; can be sterilized and used many times. Transparant body allows monitoring of dhara liquid flow and air – entrapment
3. Permanent body made of white Bronze metal with threaded end to which a disposable vasti tip can be fixed. Two sizes are available in vasti tips. Paediatric type suits all anal size.

|            |           |            |           |
|------------|-----------|------------|-----------|
| Diameter   | : 4.5 cm  | Length # 2 | : 18.5 cm |
| Length # 1 | : 19.5 cm | Length # 3 | : 19.5 cm |

## 70 MATHRA VASTI / THAILA VASTI YANTRAM



### 1. Metal Syringes

- These have Chrome plated Brass body.
- The disposable plastic vasti tip is screwed to the head.
- This tip has ideal hole size of 4 mm for ayurvedic use.
- Available in 100ml / 200 ml. sizes for Mathra and Sneha vasti procedures.
- Disposable tips come in two sizes-Adult & paediatric.

### 2. Disposable Syringes

- These are disposable vasti syringes made of high quality plastic.
- Complete suction and backward arrest even after continuous use is ensured.

# Vasti Yantra

# Basti Yantra

## 71 NIRUHA BASTI YANTRA



- Disposable tips (nethram) & vasti bags with a permanent barrel body.
- The barrel end is threaded to take up vasti bags (plastic) without leakage, using a plastic nut, through which the neck of the vasti bag is pulled out, slid over the barrel thread and nut tightened.
- 1 barrel (body) & 50 sets of disposable tips & bags of various size (500ml, 1 litre, 1.5 litre) are included in kit box.
- After 50 vasti kriyas are done with this kit, new sets of disposable tips & bags can be ordered while the body part is used permanently.
- Note: These tips can be used with Mathra & Sneha vasti yanthrams also.

## 72 UTTAR BASTI YANTRA



- Stainless steel apparatus of 10 inch length for Uttar Basti in Females.
- Adjustable to fit all sizes.
- Rubber plug provided to mark and hold the inserted position with the help of a screw.
- Uttar Basti table in stainless Steel / powder coated steel can also be provided for accurate and hygeinic treatment.

## 73 BAHYA BASTI YANTRA



GREEVA VASTHI

KATI VASTHI (B)

KATI VASTHI (S)



JAAU VASTHI (B)

JAAU VASTHI (S)

URO VASTHI

AKSHI THARPANA YANTHRA

- These are essential for executing oleation procedures like Kati Vasti, Jaanu Vasti, Greeva Vasti, Uro Vasti etc. in a fast, reliable & economical way compared to the usual Gram flour paste method or Acrylic vasti rings.
- Edges of Stainless steel Yanthra have soft rubber beading. Recyclable molding clay (included) is used to fill gaps, with fingers. Oil is poured in immediately which is not possible in any other method using Besan flour.
- Made with flexible stainless steel, these can be shaped by hand to suit the size of the body part.
- Yantra is fixed on the body by the Vulcro straps
- A Nylon Goggle with head string is supplied with the Kit as Akshi Tharpana Yanthra, for eye treatments.
- Economy set (4 piece) : Kati, Jaanu, Greeva vasti of large size plus Akshi Tharpan yanthra ( Netra Basti )
- Standard set (8 piece) : One each of the above three in large & medium sizes, urovasti, netra basti
- Premium set (11 pieces): Small, Medium & Large sizes of kati, jaanu, greeva vasti plus Urovasti & Akshi Tharpan yantra



# Other Therapy Accessories

## 74 PANCHAKARMA ACCESSORIES

Export Set



- Most of these are made of white Bronze metal, hand casted in traditional shapes and designs. Available in 3 sets. (Sets of 18 ,19 or 20 items).
- The common items in these 3 sets are dhara-Kindi with handle, Thannikinnam, Oil bowl, Bowl with handle, Thavi, Mortar & Pestle, Lotta, Polished Varpu(5"), Nasya Yanthram, Incense burner, Nilavilakku (lamp) 12", Chattukam, Polished Uruli (3"), Unpolished Uruli(12" dia), etc;
- **Econ.set: (18 items)**; with porcelain Mortar, SS nasya yanthram, Jug & Dhooma Nethram
- **Std set: (19 items)** has Shalaka, Dhooma Nethram, Mathra vasti yanthram and Niruha vasti kit
- **Export set: (20 items)** Has these additional items: Polished warp (9"), bowl, Brass spoon, Vanai and Nasya yanthram (big).

## 75 SHALAKA ( PANCHALOKHA )



- This is an indispensable tool for Agnikarma procedure.
- Made from 5 metals viz. Gold, Silver, Copper, Iron and Tin by the traditional hand casting method.
- Double tip model [Model 3] (bindhu tip & rekha tip) is used for circle & swastik procedures.
- Supplied with wooden handle.

|          |                  |
|----------|------------------|
| Length : | 28 cm / 11 inch  |
| Tip :    | 7.50 cm / 3 inch |

## 76 SHIROVASTI CAP

Standard Leather



- Adjustable leather cap made out of calf leather.
- Vinyl lining on the inner surface protects the leather from getting dirty & allows easy cleaning.
- Also available in synthetic leather in economical version .
- Size can be easily adjusted with Velcro strappings.

|        |         |
|--------|---------|
| cm :   | 80 x 20 |
| inch : | 32 x 8  |

# Other Therapy Accessories

## 77 FACIAL STEAMER / VAPORIZER



- Extra large can for long hours of continuous steam.
- Special attachments for quick relief from common cold cough & throat irritation and for soft skin.
- Large facial attachment with smooth edges for facial sauna.
- Smooth edges of nose & mouth attachment for general steam inhalation.
- Specially designed attachment for focusing steam useful for kids & elderly. Auto low water cut off.
- Quick & constant steam.

## 78 IR HEATER WITH FLAT URULI



- The 'IR' Heater is a shockproof heating system which heats oil by energy transfer with minimum electrical power.
- Since it requires a flat bottom vessel, a special Uruli of hand casted White Bronze metal is supplied with this.
- The thick metal bottom of uruli allows uniform heating of medicinal oils which usually get a burnt layer when heated in normal vessels.
- Silent,hygienic and economical heating achieved.

# Spa Furniture

## 79 SPA MASSAGE TABLE (KAYA / MAYA)



- Premium quality solid wood professional spa table.
- 2 storage drawers are provided any one of which can be positioned to access from either side of the massage table during treatment.
- The top cushion is made of 75 mm , 40 density premium grade foam upholstered with high quality leatherette that is durable & easy to clean.
- The backrest can be inclined at different positions & the reclining mechanism is made from corrosion resistant stainless steel.
- Provided with either a facial ring cushion or a breathing hole ( a face hole cutting. )
- A space can be created towards the head side to keep the decorative brass urli underneath the hole in which fragrant flowers are placed.

cm : 185 x 78 x 76  
inch : 73 x 31 x 30

## 80 MASSAGE TREATMENT TABLE



- This massage table is used for Physio & massage therapy and also as an examination table.
- Frame made of powder coated steel
- Load capacity is 200 kg.
- Bed is 75 mm thick, made of superior foam.

cm : 182x 61 x 76  
inch : 72 x 24 x 30

# Spa Furniture

## 81 KANSA VATKI BOWL & WAND



### Kansa Vataki Bowl

- The heart of the kansa vataki ayurvedic massage technique is this metal bowl.
- The bowl is specially made to help balance three doshas of human body.
- The bowl is made of 3 therapeutic metals :  
Copper: relieves pain, inflammation, arthritis  
Tin: relieves headaches and insomnia.  
Zinc: boosts immunity, digestion, good for diabetes.
- Ghee is applied on the back of the sole & other areas of the feet. Then the bowl is moved with vigorous circular & stroking movements

### Kansa Vataki Wand

- This has a hand made wooden handle attached to the kansa bowl to facilitate effortless foot massage.

## 82 FOOT BATH BOWL



- Foot bath bowls made of Copper or Brass metals.
- Wooden enclosure with wheels can also be supplied in aesthetic designs.

|           |           |
|-----------|-----------|
| cm : 35   | cm : 40   |
| inch : 14 | inch : 16 |

## 83 WOODEN MULTI PURPOSE BENCH & OIL VESSEL



- The vessel is made from single piece medicinal woods like Neem, Venga etc. by hand carving.
- This is used to collect and store Shirodhara / Pizhichil oils.
- The aesthetic wooden bench is used for placing the vessel under the Shirodhara bowl of Dharapathi as well as for seating the patient for head oiling.

## 84 WOODEN TROLLEY



- Trolley with 2 shelves and 2 drawers for therapy accessories.
- Made with seasoned wood.
- Top tray has raised edges to keep things without falling.
- Ergonomic handle for effortless movement.

# Medicinal Preparations

## 85 AYURVEDIC MEDICINES



- We also supply and export authentic ayurvedic medicines and Special, purpose made Kerala Panchakarma preparations made by reputed Kerala manufacturers as well as those prepared in the traditional way using Bronze vessels and Firewood.

Detailed Price list available on request.

## 86 AYURVEDIC OILS



- We supply and export authentic ayurvedic oils for shirodhara, abhyangam , pichu , kayasekam , nasyam , shirovasti , kativasti / janu vasti / greeva vasti etc.
- Dhara special purpose oils for relief from pain, spasm, burning sensation etc. are also supplied.
- Our shirodhara oil can be used for all persons as it is not dosha specific. Oils made by modern and traditional methods are available.

Detailed Price list available on request.

## 87 HERBAL PREPARATIONS



- We supply and export herbal preparations for all Panchakarma and Kerala specialty therapies like Potala , Swedana , Udwarthana , Lepana , Brumhana, Kshalana, Sirolepa , Shashtika , Upanaha etc.

Detailed Price list available on request.